



2015-2016 PRODUCT CATALOG

www.miulab.com

Professional manufacturer of laboratory instruments





Hangzhou Miu Instruments Co., Ltd.



MIX-1500 Micro Plate Shaker

MIX-1500 Micro Plate Shaker is mainly applied in solution mixing and shaking with PCR plate (96-holes/384-holes), cell culture plate (24-holes/48-holes/96-holes etc.) for immunoassay and dyeing experiments. It is safe used in low temperature environment or low temperature incubator. MIX-1500 is designed compact with simple and convenient operation.

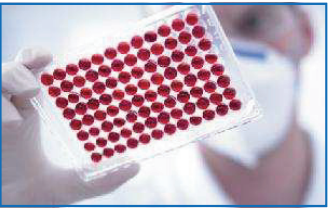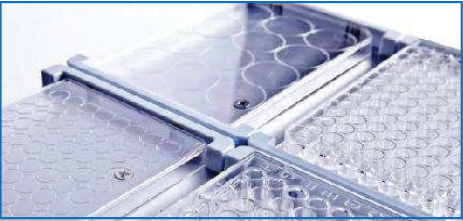

【 Product Features 】

1. LED display. Simple and easy operation designed panel. Display time and speed accurately.
2. Microprocessor controls speed and time, safe, stable, no noise.
3. Accommodate 4 standard PCR plates or micro plates, mixing trace amount of sample efficiently.
4. DC brushless motor, long life and free of maintenance. Provide both gentle and powerful shaking.
5. Timing function. Alarm after program finished.

【 Technical Parameters 】

Model	MIX-1500
Speed Range	200 ~1500rpm - standard, 200 ~1000 rpm – deep-well plate
Shaking Orbit	3mm (horizontal)
Time Range	1min ~99h59min
Sample Capacity	4 standard plates
Voltage	AC 100V ~230V, 50/60Hz
Power	50W
Fuse	250V, 1A, Φ5x20
Dimension (WxDxH)	284 x 264 x 121mm
Net Weight	5.2kgs