


JY-SPAT Horizontal Electrophoresis Tank

- Abundant buffer not only guarantees cooling effects, but also keeps the pH value stable during the whole experiment process.
- Removable electrodes enable maintenance and replacement easy and quick.
- With background color design on the tray, convenient to adding samples.
- It's more convenient to open the lid with the easy-open buttons.
- It prevents buffer volatilizing and electricity leakage with transparent lid.
- Auto-switch-off when the lid is moved.

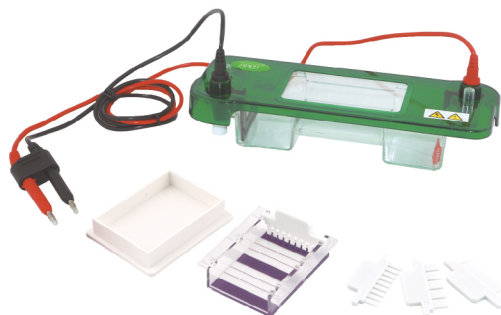

JY-SPAT Packing List


Parts name	Quantity
Lid with cables	1
Body tank with electrodes	1
48×75mm Gel tray	1
Gel box	1
1.0mm thickness 3-well comb	2
1.5mm thickness 3-well comb	2
1.0mm thickness 5-well comb	2
1.5mm thickness 5-well comb	2
1.0mm thickness 9-well comb	2
1.5mm thickness 9-well comb	2

JY-SPAT Technical Specifications

Gel size (W×L)	48×75 (mm)
No. of samples	3~18
Total buffer volume	150 (ml)
Dimension (L×W×H)	250×90×50 (mm)
Weight	0.4 (kg)

JY-SPBT Horizontal Electrophoresis tank

- Abundant buffer not only guarantees cooling effects, but also keeps the pH value stable during the whole experiment process.
- Removable electrodes enable maintenance and replacement easy and quick.
- With background color design on the tray, convenient to adding samples.
- It's more convenient to open the lid with the easy-open buttons.
- It prevents buffer volatilizing and electricity leakage with transparent lid.
- Auto-switch-off when the lid is moved.

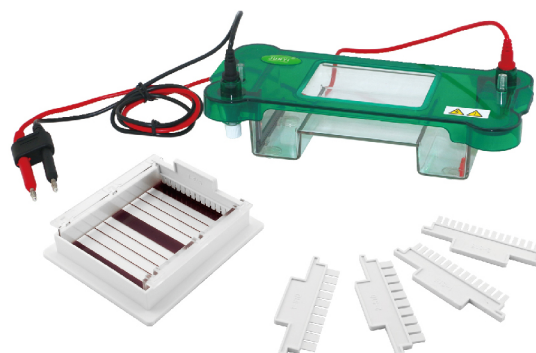

JY-SPBT Packing List

Parts name	Quantity
Lid with cables	1
Body tank with electrodes	1
78×100mm Gel tray	1
Gel box	1
1.0mm thickness 10-well comb	2
2.0mm thickness 10-well comb	2
1.0mm thickness 15-well comb	2
2.0mm thickness 15-well comb	2

JY-SPBT Technical Specifications

Gel size (W×L)	78×100 (mm)
No. of samples	10~30
Total buffer volume	400 (ml)
Dimension (L×W×H)	365×105×55 (mm)
Weight	0.6 (kg)