


JY-SPFT Horizontal Electrophoresis Tank

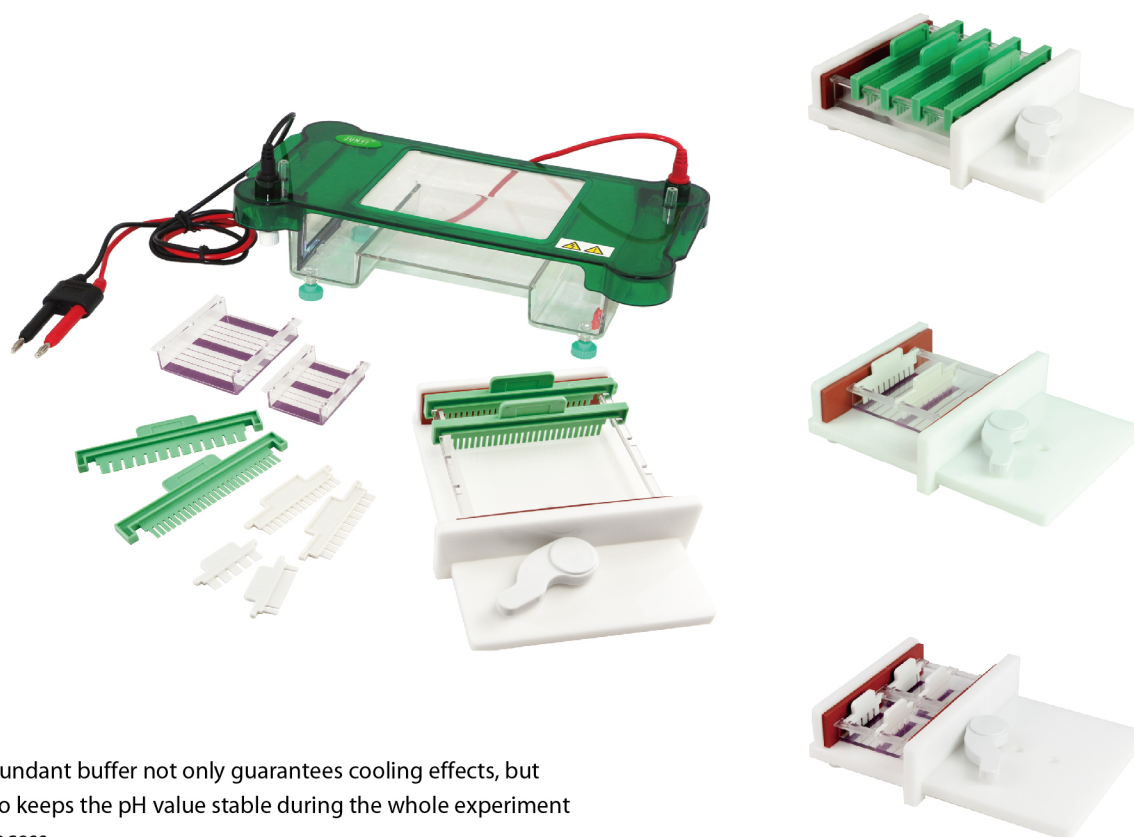

- Abundant buffer not only guarantees cooling effects, but also keeps the pH value stable during the whole experiment process.
- Removable electrodes enable maintenance and replacement easy and quick.
- With the fluorescent ruler on the tray, can be observed directly into the ultraviolet.
- With a special gel casting, preventing the gel leakage, operation is simple and convenient.
- It's more convenient to open the lid with the easy-open buttons.
- It prevents buffer volatilizing and electricity leakage with transparent lid.
- Auto-switch-off when the lid is moved.

JY-SPFT Technical Specifications

Gel size (W×L)	130×130 (mm)
	78×100 (mm)
	48×75 (mm)
No. of samples	3~104
Total buffer volume	600 (ml)
Dimension (L×W×H)	305×152×80 (mm)
Weight	1.7 (kg)

JY-SPFT Packing List

Parts name	Quantity
Lid with cables	1
Body tank with electrodes	1
130×130mm Gel tray	1
78×100mm Gel tray	1
48×75mm Gel tray	2
Gel casting stand	1
1.0mm thickness 13-well comb	2
1.5mm thickness 13-well comb	2
1.0mm thickness 26-well comb	4
1.5mm thickness 26-well comb	2
1.0mm thickness 10-well comb	2
1.0mm thickness 15-well comb	2
1.5mm thickness 3-well comb	2
1.0mm thickness 5-well comb	2
level adjusting button	4