


## JY-Mini-P4 Vertical Electrophoresis Tank



- With two gel castings stands, 1-4 plates of gels can be made at the same time.
- The installation of glass plates can be completed within 15s, easy and quick.
- With a gel casting stand in original position, no need to move glasses from gel making to electrophoresis running. It is more convenient to observe gel preparing through both sides of the glasses.
- With the background color, it is easy to observe sample loading and electrophoresis running.
- Abundant buffer not only guarantees a cooling environment, but also keeps the pH value stable during the whole experiment.
- With protective strip of positive electrode, it can prevent producing bubbles and also protect platinum silk from damage.
- Auto-switch-off when the lid is moved.

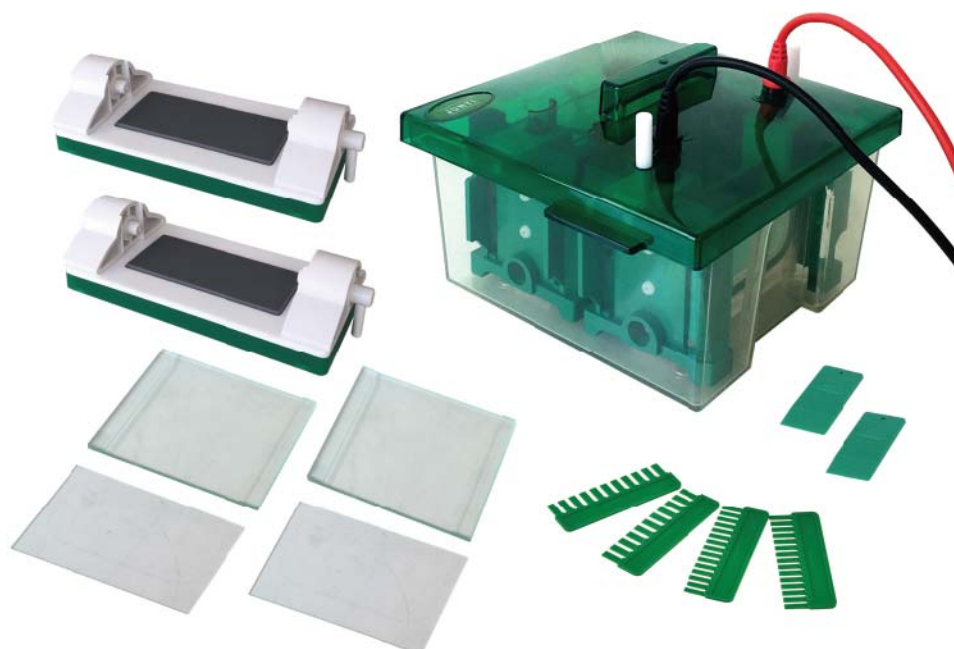

## JY-Mini-P4 Technical Specifications

Glass plate size (W×L)	101×82 (mm)
Gel size (W×L)	101×73 (mm)
Gel thickness	1.0 (mm)
Number of gel	1~4 pc(s)
No. of samples	11~60
Total buffer volume	1200 (ml)
Dimension (L×W×H)	175×175×100 (mm)
Weight	2.5 (kg)

## JY-Mini-P4 Packing List

Parts name	Quantity
Lid with cables	1
Body tank	1
Gel casting core with electrodes	2
Gel casting stand	2
Short glass plates	6
Spacer plate with 1.0mm glass spacer	6
1.0mm thickness 11-well comb	6
1.0mm thickness 15-well comb	6
Gel releaser	2