

# MULLAB

## Mini Centrifuge

Fully transparent round upper cover, equipped with multiple rotors to meet various experimental needs



Mini-10K+C

# Centrifugal force increased by 1.7 times

Large radius rotor with 8 tube holes, centrifugal force is 1.7 times bigger than ordinary 6-hole rotor at the same rolling speed.

Fully transparent cover

32-hole rotor

LED display



Mini-4KC/6KC/7KC/10KC

# The upper cover can automatically pop open 95°

The cover press-button is outward designed with start/stop function. Press the button, the cover will open automatically to 95°, which is very convenient for one-hand operation.

## Upper cover switch





Mini-10K+C



# Unique rotor buckle design

Unique rotor buckle design, convenient for rotor replacement, take and put.

Mini-4KC/6KC/7KC/10KC/10K+C

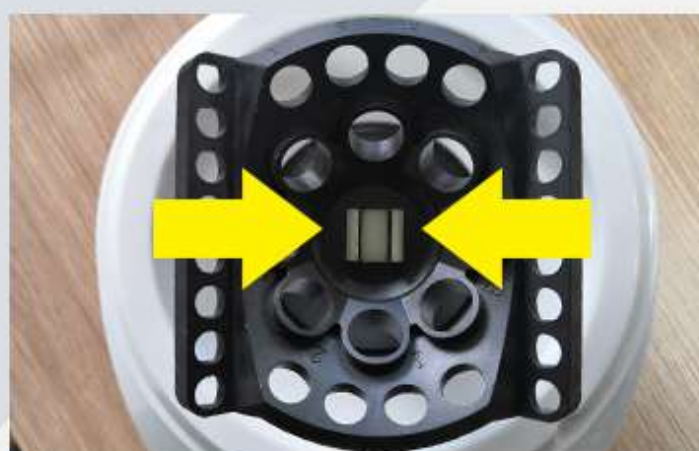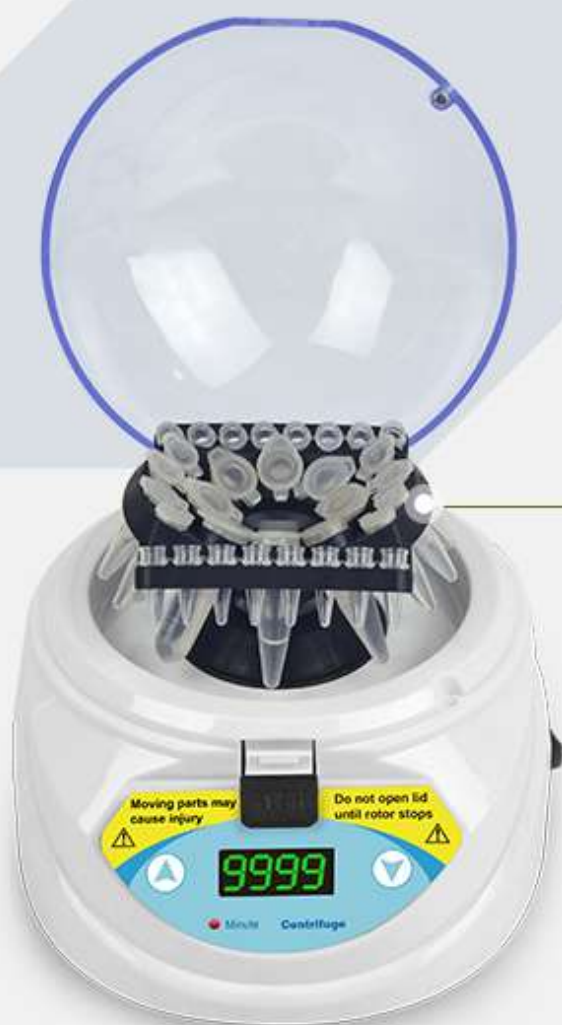

Squeeze the slot from both sides to the middle to take out the rotor and replace the rotor, which is convenient and quick.

Combined rotor & tube adaptor



Combined rotor



Tube adaptor

# Ultra-low noise, durable motor

Working is silent and stable. Free-maintenance motor, long-life and safe.

Mini-4KC/6KC/7KC/10KC/10K+C



Mini-10K+C

# Technical Parameter

Mini-4KC/6KC/7KC/10KC/10K+C



Model	Mini-4KC	Mini-6KC	Mini-7KC	Mini-10KC	Mini-10K+C
Speed Range	4000rpm	6000rpm	7000rpm	10000rpm	3krpm ~10krpm (increment:1000rpm)
RCF	approx.1000g	approx.2300g	approx.3100g	approx.5600g	approx.500g ~5600g
Tube adapter	8x2.0ml/1.5ml centrifuge tubes, 8x0.5ml and 8x0.2ml, 2 PCR strips or 2x8x0.2ml tubes	8x2.0ml/1.5ml centrifuge tubes, 8x0.5ml and 8x0.2ml, 2 PCR strips or 2x8x0.2ml tubes	8x2.0ml/1.5ml centrifuge tubes, 8x0.5ml and 8x0.2ml	8x2.0ml/1.5ml centrifuge tubes, 8x0.5ml and 8x0.2ml	8x2.0ml/1.5ml centrifuge tubes, 8x0.5ml and 8x0.2ml, 2 PCR strips or 2x8x0.2ml tubes (speed<=6000rpm)
Time Range	1s ~9999min or continuous	1s ~9999min or continuous	1s ~9999min or continuous	1s ~9999min or continuous	1s ~999min or continuous
Working Noise	≤45dB	≤45dB	≤55dB	≤55dB	≤55dB
Voltage	AC 220V / AC 110V, 50/60Hz	AC 220V / AC 110V, 50/60Hz	AC 220V / AC 110V, 50/60Hz	AC 220V / AC 110V, 50/60Hz	AC 220V / AC 110V, 50/60Hz
Power	≤45W	≤45W	≤45W	≤45W	≤45W
Fuse	250V, 1A, Φ5x20	250V, 1A, Φ5x20	250V, 1A, Φ5x20	250V, 1A, Φ5x20	250V, 1A, Φ5x20
Dimension (WxDxH)	W.156 x D.176 x H.121mm	W.156 x D.176 x H.121mm	W.156 x D.176 x H.121mm	W.156 x D.176 x H.121mm	W.156 x D.176 x H.121mm
Net Weight	1.5kgs	1.5kgs	1.5kgs	1.5kgs	1.5kgs

HANGZHOU MIU INSTRUMENTS CO.,LTD.

PROFESSIONAL MANUFACTURER OF LABORATORY INSTRUMENTS